Experiment 2 Verification of Newton's Second Law of Motion by Atwood Machine.

Physics Lab 1
Spring 2021-22
Department of Physics
American International University-Bangladesh

Objectives:

To establish the relationship between force and acceleration, thus verify Newton's second law of motion.

Outcomes:


After completing this experiment student should be able to answer the following questions:

- What is the relationship between force and acceleration for an object according to Newton's second law of motion?
- What is the basic concept of net force?
- How an Atwood machine can be constructed? How different forces and acceleration work for the Atwood machine.
- Why the experimental accelerations vary from the theoretical accelerations?
- What is the meaning of a linear relationship and how it looks in a graph?

Theory: Newton's Second Law of Motion


Newton's second law of motion says that FORCE equals MASS times ACCELERATION



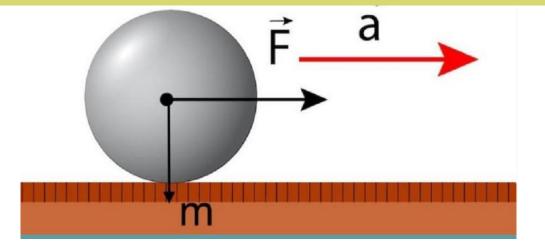

For a particular mass:

 $Acceleration \propto Force$

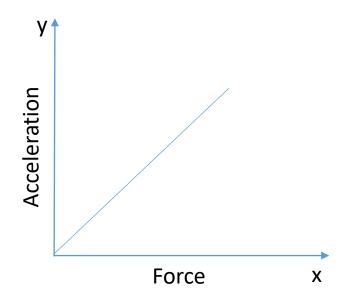

Atwood Machine: Theory

Applying Newton's 2^{nd} Law: $F_{net} = ma$

• For M:

$$\mathbf{F}_{\text{net}} = \mathbf{T} - \mathbf{M}\mathbf{g} = -\mathbf{M}\mathbf{a},$$

• For m:


$$\mathbf{F}_{\text{net}} = \mathbf{T} - \mathbf{mg} = \mathbf{ma}$$

• Solving these two equations: the theoretical acceleration

$$a_{th} = \frac{g}{(M+m)} (M-m)$$

• Keeping (**M**+**m**) constant at any particular place (**g** is constant), we get


$$a_{th} \propto (M-m)$$

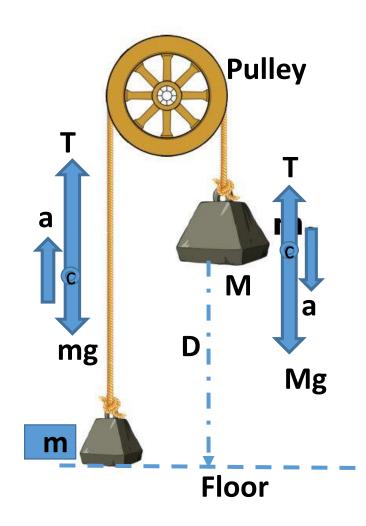

Atwood Machine: Experiment

- M falls a distance D in time t from rest.
- Applying the knowledge of equations of motion (D=ut+1/2at²), we get the experimental acceleration

$$a_{ex} = \frac{2D}{t^2}$$

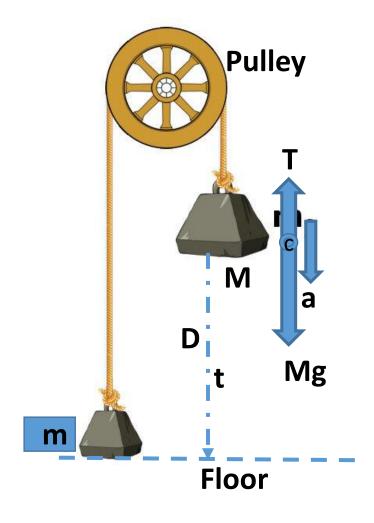

Verification of Newton's 2nd Law

- Atwood Machine: Theory Predicts

$$a_{th} \propto (M-m)$$

Atwood Machine: Experimental Result


$$\mathbf{a}_{\mathrm{ex}} = \frac{2\mathbf{D}}{\mathbf{t}^2}$$

If we find

$$a_{ex} \propto (M - m)$$

we can say, Newton's 2nd law is verified.

 $Acceleration \propto Mass \ difference$

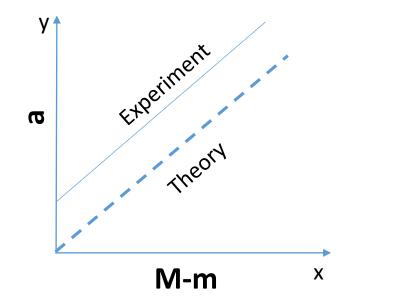

Apparatus


Atwood Machine:

Pulley, two hangers, different masses, string, stand and clamp.

Measurement of D and t:

Meter scale and stop watch.



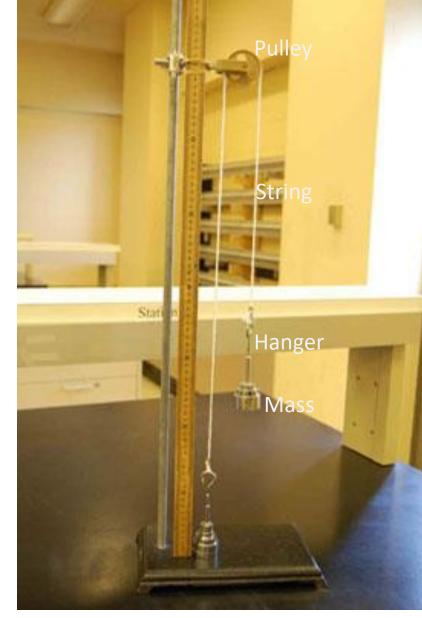

Procedure:

- Hold the lighter mass on the floor attached to one end of a string. The
 heavier one attached to the other end of the string is up in the air at a
 height D from the floor. Measure D with a meter scale.
- Now release the lighter mass and measure the time the heavier mass takes to fall onto the floor. Run the experiment for 7 different mass-differences, (M m). For each run, obtain the value of the acceleration in (m/s^2) experimentally as well as theoretically. Make sure to keep total mass (M + m) always constant.
- Using Excel plot acceleration (a_{th} and a_{ex}) versus mass difference (M m) graph.

Video lecture on Procedure:



Lab Works:

• Complete the data table with the calculations.

• Draw the acceleration vs mass difference graph in Excel. Plot both the accelerations (theoretical and experimental) on the same graph paper.

Analyze the result.

Discussion on Outcomes of the Lab

- What is the relationship between force and acceleration for an object according to Newton's second law of motion?
- What is the basic concept of net force?
- How an Atwood machine can be constructed? How different forces and acceleration work for the Atwood machine.
- Why the experimental accelerations vary from the theoretical accelerations?
- What is the meaning of a linear relationship and how it looks in a graph?